


# Positioning of Primus Dryers



## Single Installation

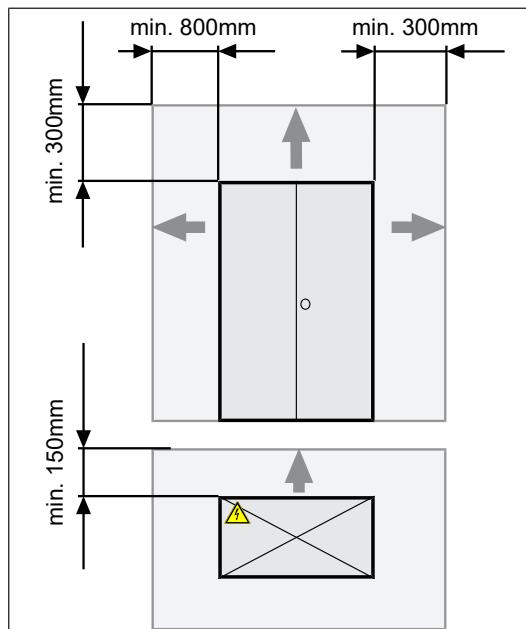

The installation site must be sufficiently level. Ensure that the system stands firmly in a 100 % vertical position and does not wobble.

- Ensure that the doors align perfectly at the front. If necessary, suitable leveling material must be placed under the drying locker if the floor is uneven.




**CAUTION** During mounting, ensure that there is unobstructed ventilation at the rear of the device and access to the main switch and to the mains connection at the rear left on the top of the drying system.

Observe the following minimum distances from the walls and ceiling.

| Space requirement/minimum distance |        |        |        |
|------------------------------------|--------|--------|--------|
| Ceiling/above                      | Left   | Right  | Rear   |
| 300 mm                             | 800 mm | 300 mm | 150 mm |

## Installation in a Row

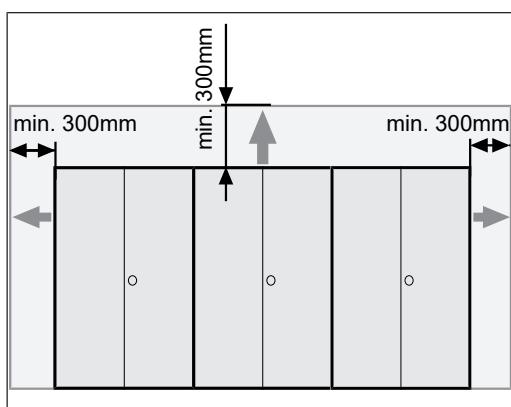


### Spaces when several drying lockers are installed in a row

If several drying lockers are installed next to each other in a row, the distance between the individual drying lockers can be reduced.



**NOTE** There must be at least 300 mm from adjoining walls and also from the ceiling, and at least 150 mm from the back wall.

## 5.3.3 Primus Set 2 Wall Holder

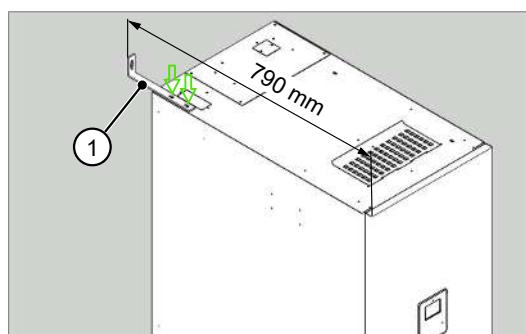

**CAUTION** Free-standing drying lockers could fall over and injure persons.

1. Secure the wall holder [1] to the top of the drying locker using 2 screws.
2. Mount the drying locker on the wall using a suitable means of fastening.



**NOTE** The Primus Set 2 does not require 150mm clearance at the rear, this is set by the Wall Holder.